

ACTIVHEAL® DISCOVER MORE.



DISCOVER MORE ABOUT ACTIVHEAL®

ACTIVHEAL® WAS ORIGINALLY DEVELOPED BY ADVANCED MEDICAL SOLUTIONS TO OFFER THE NHS A MORE AFFORDABLE, CLINICALLY EFFECTIVE WOUND CARE DRESSING RANGE.

The ActivHeal® range was initially developed in 2005, specifically to meet the needs of the NHS, by offering cost savings without compromising on clinical outcomes or patient care. We continue to offer clinicians a simple, clear and cost conscious wound care range whilst ensuring the best clinical outcome for the patient. ActivHeal® offers a full range of wound care dressings. The ActivHeal® dressing range facilitates moist wound healing, ensuring that all patients' wounds have the optimal environment to progress. Available in a range of sizes and shapes, each dressing type offers a variety of features to manage both chronic and acute wounds.

BUT ACTIVHEAL DOESN'T STOP THERE.

- ▶ Over 30 years of UK manufacturing
- ▶ Clinically proven performance at prices you may find surprising
- ▶ Fully supported with Society of Tissue Viability endorsed higher tier education
- ▶ Simplified packaging helping you to get it right first time, every time



EDUCATION SUPPORT

First to offer STV endorsed Higher Tier Education platform



SIMPLIFIED PACKAGING

Getting it right first time, every time



HIGH QUALITY

Clinically proven performance



MORE AFFORDABLE

Clear and transparent pricing



MADE IN THE UK

Over 30yrs of UK Manufacturing





Find out more at
academy.activheal.com





DISCOVER
MORE

ABOUT
ACTIVHEAL®

4



ACTIVHEAL® PHMB FOAM
SILICONE BORDER

10



ACTIVHEAL®
FOAM PHMB

12



ACTIVHEAL®
AQUAFIBER™ EXTRA

14




ACTIVHEAL®
AQUAFIBER™ AG

16



ACTIVHEAL®
SILICONE FOAM

18

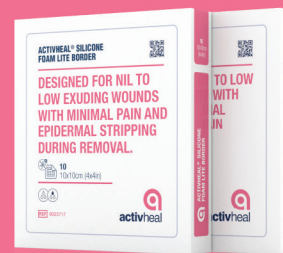

ACTIVHEAL®
SILICONE FOAM LITE

19



ACTIVHEAL®
AQUAFIBER™

20



ACTIVHEAL®
ALGINATE

21



ACTIVHEAL®
FOAM RANGE

22



ACTIVHEAL® SILICONE
WOUND CONTACT LAYER

24



ACTIVHEAL®
HYDROGEL

25



DISCOVER
MORE

NOW THAT'S VERY CLEVER.

HIGH-QUALITY, MULTI-TIERED EDUCATION
TO SUPPORT YOU IN YOUR ROLE



ACCESS OUR MULTI-TIERED WOUND CARE EDUCATION PROGRAMME.

WHAT IS THE ACTIVHEAL® ACADEMY?

The ActivHeal® Academy was developed to provide Health Care Professionals with accessible, non product related wound care education.

The aim is to assist Health Care Professionals with development of their wound care knowledge, that will lead to achieving positive outcomes for patients.

The ActivHeal® Academy is an online learning platform, with nine simplified guides and five higher-tier modules. The courses are freely available for Health Care Professionals to assist in their day to day practice while supporting their Continuing Professional Development.

EDUCATION SUPPORT

Working in partnership with wound care organisations to ensure that all ActivHeal® Academy users have access to a higher level of quality assured education that has been reviewed by specialists in wound care. The ActivHeal® Academy higher modules were the first to receive their endorsement, (The Society of Tissue Viability - STV), ensuring that our education provides:

High-quality education that meets STV endorsement standards.

All material meets stated learning objectives.

Educational material that is free from bias or product promotion.

Access to the general wound care community.



96%

**OF CLINICIANS SAID THE PROVISION
OF EDUCATION MATERIALS IS
IMPORTANT TO THEM***

NINE SIMPLIFIED GUIDES

Covering a variety of wound care topics from the process of wound healing to pressure ulcer management and are available to view online to download as a PDF, printed form and presentation format.

FIVE COMPREHENSIVE HIGHER MODULES

Based on a higher learning university level and currently endorsed by the Society of Tissue Viability and South African Medical Association

7TH EDITION EDUCATION POCKET GUIDE TO WOUND CARE

164 page wound care education Guide. Packed with useful information for Health Care Professionals and is exclusively available to ActivHeal® customers.

- ✓ Online tests and PDF certificates
- ✓ Supports Continuing Professional Development (CPD)
- ✓ Facility to create groups
- ✓ Optimised for use on mobiles and tablets



To discover more about the ActivHeal® Academy Society of Tissue Viability endorsed higher modules and to access our full range of nine simplified guides, register online at:

academy.activheal.com

SIGN UP TODAY

THIS PRACTICAL WOUND CARE GUIDE IS PACKED FULL OF USEFUL INFORMATION, IT'S ALSO SMALL ENOUGH TO FIT INTO YOUR POCKET OR DIARY

7TH EDITION

activheal ACADEMY

MODULE 1: ANATOMY & PHYSIOLOGY OF SKIN

MODULE 2: WOUND HEALING

MODULE 3: PATIENT AND WOUND ASSESSMENT


MODULE 4: WOUND CLASSIFICATION AND WOUND MANAGEMENT



DISCOVER
MORE

OH YEAH, OBVIOUSLY!

EASY TO USE, COLOUR-CODED
PACKAGING TO SAVE YOU TIME

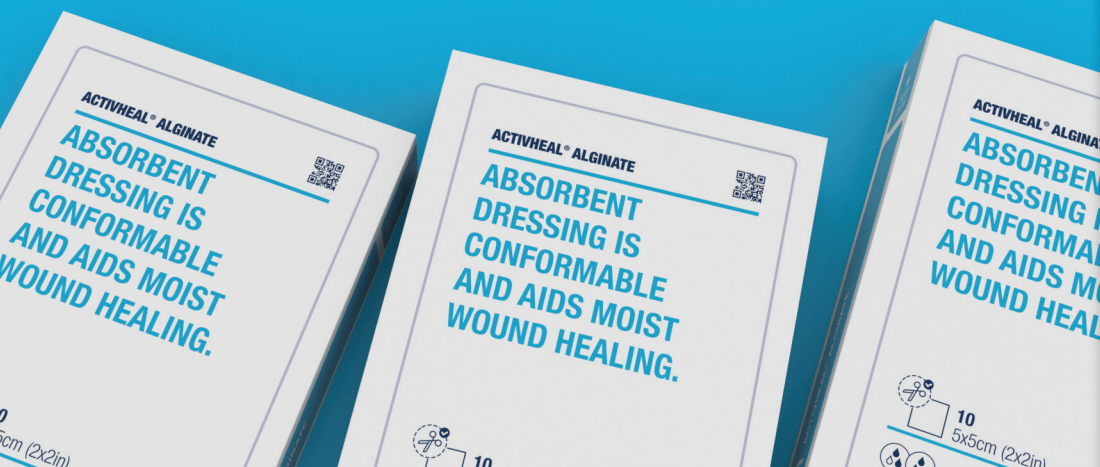

SIMPLIFIED PACKAGING GETTING IT RIGHT FIRST TIME, EVERY TIME.

We understand that choosing the right dressing for the right wound type can save Health Care Professionals' time and valuable resources. All of our products are innovatively designed and packaged to aid in selecting dressings, allowing Health Care Professionals to spend more time with their patients as well as improving clinical outcomes.

WHAT CHANGES HAVE BEEN MADE TO THE ACTIVHEAL® PACKAGING TO HELP YOU?

Our simplified approach to packaging incorporates colour coded boxes with generic names and clear descriptions to help you make accurate selections and assist in reducing unnecessary waste.



Improved inner packaging and
Instructions for Use (IFU).



QR code and messaging for extra support.



Easy to identify icons for absorbency
and cutting usage.



New simplified dressing shape icons for
everyday use.

GREAT SAVINGS WITHOUT COMPROMISING ON PATIENT CARE.

ActivHeal® wound care includes every day and advanced dressings ensuring your wound care needs are met in one range. But that's not all. Our balanced, clear and transparent pricing can provide savings across both acute and community health care settings. ActivHeal® wound care range is clinically effective, cost conscious, meaning you can save valuable resources with no impact on patient outcomes.

SWITCH TO ACTIVHEAL® FOR CLINICALLY PROVEN PERFORMANCE AT PRICES YOU MIGHT FIND SURPRISING.



DISCOVER
MORE

THAT'S UNEXPECTED!

AFFORDABLE WOUND CARE WITHOUT
COMPROMISING ON PATIENT CARE



 **MORE AFFORDABLE**

activheal.com



DISCOVER
MORE

SO SUPER RESOURCEFUL!

FREE ACCESS TO CLINICAL RESOURCES
TO SUPPORT YOU IN USING
OUR PRODUCTS



CLINICALLY PROVEN WOUND CARE PRODUCTS, SUSTAINABLY MADE IN THE UK.

ActivHeal® products are developed and manufactured in the UK. The ActivHeal® wound care dressing range has undergone extensive research and development to support the healing process at each stage, providing the optimal conditions to aid effective healing. Our dressings are clinically proven to deliver high-performance when used correctly.

CLINICAL RESOURCES


Our clinical resources include case studies, clinical evaluations and product factsheets to help you discover more about our product performance.

98%

OF USERS RATED ACTIVHEAL®
PRODUCTS TO BE OK, GOOD VERY
GOOD OR EXCELLENT*

Our customers agree that ActivHeal® products are as good as or better than our leading competitors.*

Image taken from our PHMB Mode
of Action video. Link below.

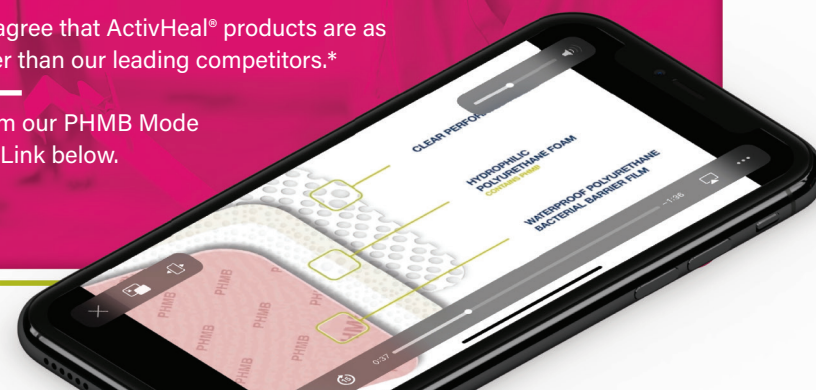

DRIVING QUALITY WHILST MAINTAINING LOW COST.

Our experienced research and development team (R&D) and ongoing investment into our manufacturing site continues to drive our high quality products and allows us to expand product ranges whilst maintaining a cost conscious approach.

- ✓ MDSAP certified
- ✓ Quality Management Certification (ISO 13485:2016)
- ✓ Onsite Quality and Regulatory team
- ✓ Dedicated in-house Research and Development department and Clinical Research team

RELIABLE SUPPLY CHAIN

Our UK production facility and central location ensures a high completion rate and on-time delivery of orders. An outstanding industry performance.

ADVANCED MEDICAL SOLUTIONS HAS COMMITTED TO:



Reduce greenhouse gas

Through effective energy management and initiatives

Creative thinking

Encourage our people to bring forward new initiatives

Sustainable sources

Ensuring our products are from sustainable sources where possible

Recyclable materials

Working with our suppliers to use recyclable materials wherever possible

Carbon footprint

Reduction in business miles by combining deliveries



**VERY
SENSIBLE!**
DRIVING QUALITY WHILE
MAINTAINING LOW COST



ACTIVHEAL® PHMB FOAM SILICONE BORDER A SOFT, CONFORMABLE DRESSING WHICH IS EFFECTIVE AGAINST GRAM NEGATIVE BACTERIA, GRAM POSITIVE BACTERIA AND YEAST.

KEY FEATURES:

// High absorbency // Soft, flexible and conformable // Supports non traumatic removal // Prevents colonisation and proliferation of bacteria // Suitable for moderate to heavily exuding wounds

KEY BENEFITS

ActivHeal® PHMB Foam Silicone Border dressing contains the antimicrobial substance polyhexamethylene biguanide (PHMB, polihexanide), which kills and inhibits the growth of bacteria. The PHMB released in the presence of wound exudate is an effective antimicrobial agent against a broad spectrum of microorganisms frequently associated with the bacterial colonisation and infection of wounds, for up to seven (7) days, based on in-vitro testing. The perforated wound contact layer contains a gentle silicone adhesive that provides secure, non-irritating adhesion and supports non-traumatic removal during dressing changes. The low friction, waterproof polyurethane film provides a bacterial barrier to the wound. The perforated wound contact layer contains a gentle silicone adhesive that provides secure, non-irritating adhesion and supports non-traumatic removal during dressing changes.

High absorbency - Prevents excess fluid causing maceration to healthy skin and can reduce the number of dressing changes

Soft, flexible and conformable - comfortable to wear and can to adjust to difficult anatomical areas

Waterproof and bacterial barrier - can shower with the dressing in situ

Perforated wound contact layer with a gentle silicone adhesive perforated film - protects the wound bed from adhering to the dressing



Secure but gentle adhesion - allows the dressing to stay in place and minimise epidermal stripping

Supports non traumatic removal - assists with minimising trauma at dressing change

Effective antimicrobial - kills and inhibits the growth of bacteria and prevents colonisation and proliferation of bacteria

Promotes a moist wound environment - assists with wound healing progress

Manages wound exudate and microbial contamination - which reduces wound pain during wound healing.

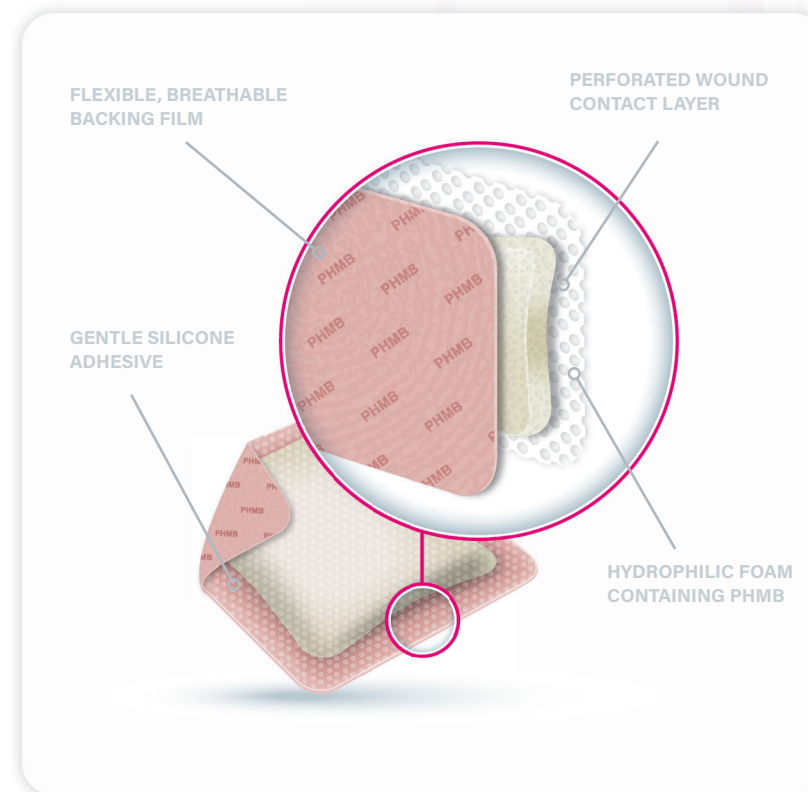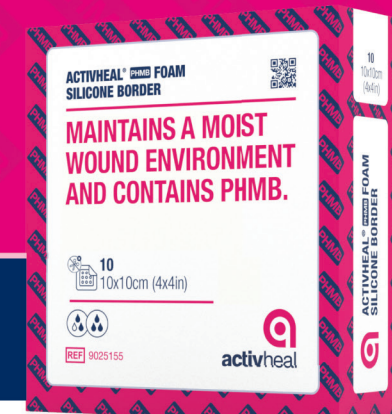

Effective Performance

ActivHeal® PHMB Silicone Foam is an effect antimicrobial dressing. PHMB is released in the presence of wound exudate and is an effective antimicrobial agent against a broad spectrum of microorganisms¹:

Challenged Organisms

- ✓ **Escherichia coli**
- ✓ **Klebsiella pneumonia**
- ✓ **Streptococcus pyogenes**
- ✓ **Candidia albicans**
- ✓ **Kelbsiella pneumonia**
- ✓ **Staphylococcus aureu**
- ✓ **Serratia marcescens**
- ✓ **MRSA**
- ✓ **VRE**
- ✓ **Streptococcus pyogenes**
- ✓ **Pseudonymous aeruginosa**
- ✓ **Serratia marcescens**



**KILLS 99.99% OF BACTERIA, YEASTS
AND MOULD WITHIN 24 HOURS²**



**BARRIER TO
BACTERIAL
PENETRATION**



**HELPS PREVENT
COLONISATION
AND PROLIFERATION**



**NON TRAUMATIC
REMOVAL**



**GENTLE SILICONE
ADHESIVE**

INDICATIONS FOR USE

ActivHeal® PHMB Foam Silicone Border dressing is indicated for moderately to heavily exuding acute or chronic wounds that are infected or at risk of infection

- ✓ **Surgical wounds**
- ✓ **Pressure ulcers**
- ✓ **Leg ulcers**
- ✓ **Diabetic ulcers**



SIZES AND CODES

ACTIVHEAL® PHMB FOAM SILICONE

SIZE (CM)	QTY	PRODUCT CODE	NHS SUPPLY CHAIN	DT PIP CODE
8x8	10	9027449	ELA1424	4198784
10x10	10	9027456	ELA1422	4198826
12.5x12.5	10	9027463	ELA1425	4198792
10x20	10	9027470	ELA1423	4198842
15x15	10	9027487	ELA1426	4198800
20x20	10	9027494	ELA1428	4198776
Sacral (18.5x19.5)	10	9027500	ELA1427	4198818

ACTIVHEAL® PHMB FOAM NON-ADHESIVE CONTAINS AN ANTIMICROBIAL WHICH KILLS AND INHIBITS THE GROWTH OF A BROAD SPECTRUM OF MICROORGANISMS.

KEY FEATURES:

// Reduces the risk of maceration // Sustained antimicrobial release and effectiveness for up to 7 days // Broad spectrum antimicrobial action on MRSA, VRE, Pseudomonas aeruginosa, Candida albicans // Helps prevent recolonisation



KEY BENEFITS

ActivHeal® Non Adhesive PHMB Foam dressings are sterile antimicrobial wound dressings, consisting of a hydrophilic polyurethane foam which is designed to absorb exudate. The foam contains the antimicrobial substance polyhexamethylene biguanide (PHMB, polihexanide), which kills and inhibits the growth of bacteria.

High absorbency of exudate - that lessens the risk of leakage and can reduce the number of dressing changes

Reduces the risk of maceration - helps to protect the wound edges and surrounding skin.

Promotes healing through a moist wound environment - assists with wound healing progression

Sustained antimicrobial release and effectiveness for up to 7 days - helps prevent recolonisation

Effective against MRSA, VRE, Pseudomonas aeruginosa, Candida albicans³ - helps to reduce local infection caused by these organisms

Conformable - adaptable to awkward areas


INDICATIONS FOR USE

Indicated for moderate to heavily exuding chronic and acute wounds that are infected or are at risk of infection. Suitable to use under compression bandaging. The dressing may be used during the healing process on:

- ✓ Leg and foot ulcers
- ✓ Pressure ulcers
- ✓ Diabetic ulcers
- ✓ Surgical wounds

HYDROPHILIC ANTIMICROBIAL PHMB
POLYURETHANE FOAM PAD


WATERPROOF
POLYURETHANE
BACTERIAL
BARRIER FILM

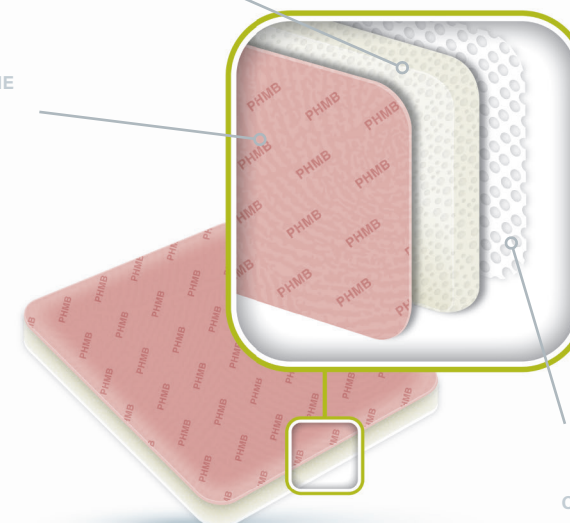

PERFORATED
WOUND
CONTACT LAYER

Challenge Organism 3

Total eradication of challenged organisms within 24 hours.

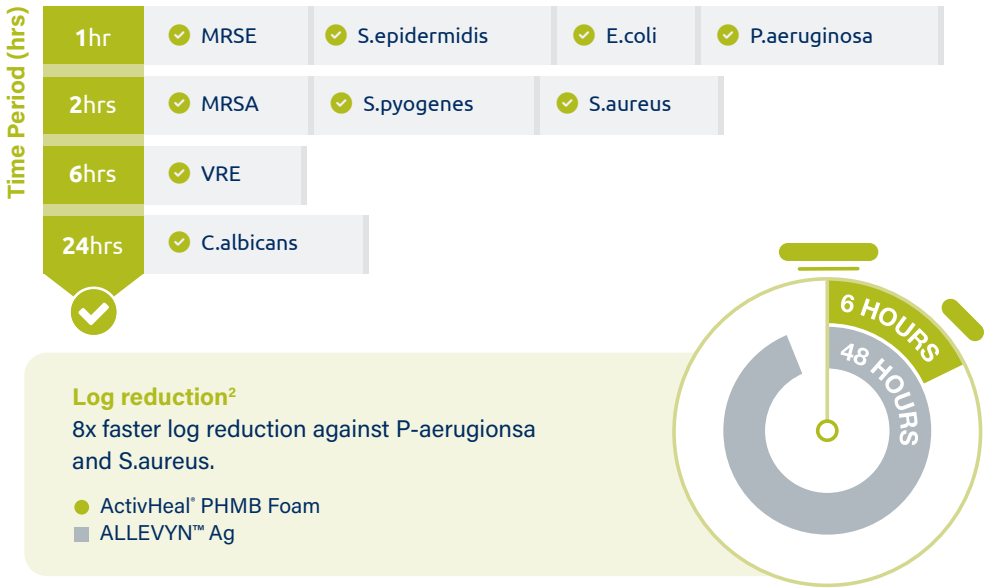

Performance

Excellent fluid handling reduces the number of dressing changes, increases wear time and reduced maceration¹.

Non-Border PHMB Foam
31g
10cm²/24hr

ALLEVYN[™] Ag
18g
10cm²/24hr

Total Fluid Handling Performance¹
Be confident this foam can handle patient exudate.

SIZES AND CODES

ACTIVHEAL [®] PHMB FOAM NON-ADHESIVE				
SIZE (CM)	QTY	PRODUCT CODE	NHS SUPPLY CHAIN	DT PIP CODE
5x5	10	9027517	ELA788	404-2024
7,5x7,5	10	9027524	ELA789	404-2032
10x10	10	9027531	ELA790	404-2040
10x20	10	9027548	ELA791	404-2057
15x15	10	9027555	ELA792	404-2065
20x20	10	9027562	ELA793	404-2073



ACTIVHEAL® AQUAFIBER™ EXTRA CREATES A MOIST WOUND HEALING ENVIRONMENT FOR MODERATE TO HEAVILY EXUDING WOUNDS.

KEY FEATURES:

// Highly absorbent // High wet tensile strength // High gelling // Reduces risk of maceration // Aids autolytic debridement




KEY BENEFITS


ActivHeal® AquaFiber™ Extra is an innovative combination of two proven wound care fiber technologies. It is a soft, conformable, highly absorbent non woven dressing, with a unique patented quilting design which turns into a soft, cohesive gel when in contact with wound exudate. This assists with autolytic debridement and maintains a moist wound healing environment following the principles of advanced wound care. ActivHeal® AquaFiber™ Extra effectively addresses the key challenges of exudate management with improved absorption, gelling and tensile performance.

INDICATIONS FOR USE

ActivHeal® AquaFiber™ Extra is indicated for the management of the following moderately to heavily exuding acute and chronic wounds:


- ✓ Pressure ulcers
- ✓ Venous leg ulcers
- ✓ Arterial leg ulcers
- ✓ Diabetic ulcers
- ✓ Graft and donor sites
- ✓ Post-operative surgical wounds
- ✓ Cavity wounds
- ✓ Superficial and partial thickness burns

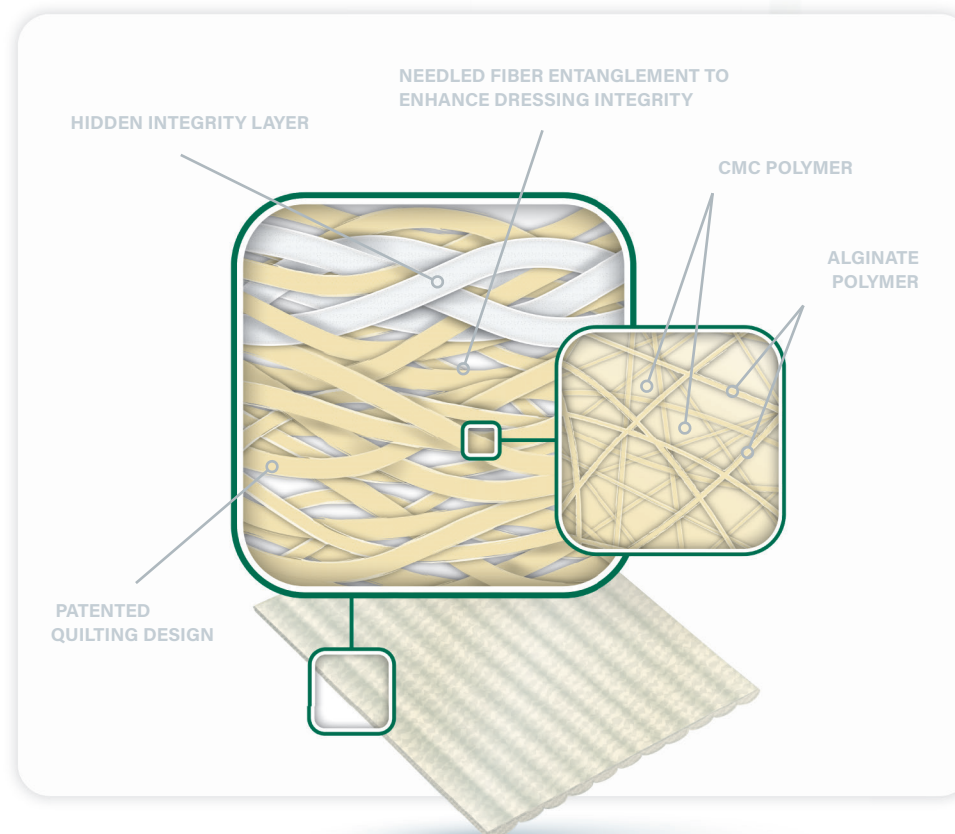

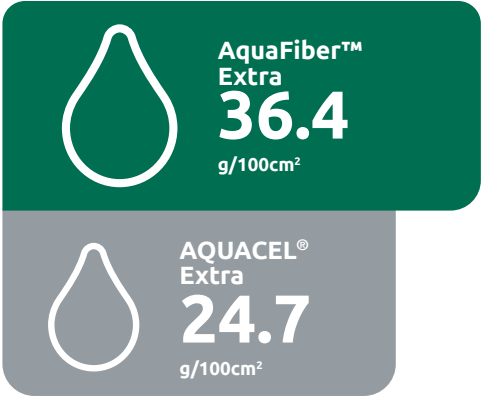

High absorbancy

The unique patented quilting design & fiber engineering of ActivHeal® AquaFiber™ Extra allows the dressing to absorb exudate & transform into a cohesive gel.

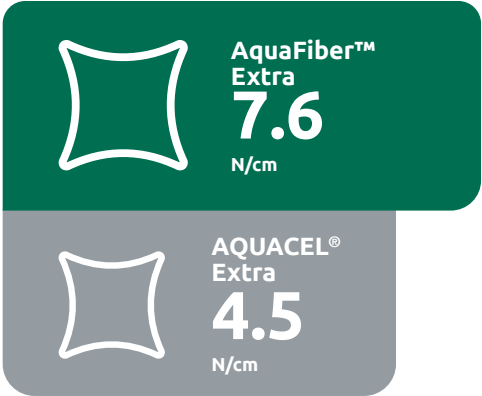

High wet tensile strength

The unique combination of the patented quilting design and hidden integrity layer ensures intact removal of the ActivHeal® AquaFiber™ Extra dressing.

Increased vertical wicking
In-house analysis¹ has confirmed that the vertical wicking of the ActivHeal® AquaFiber™ Extra can be potentially helpful in the prevention of maceration of the peri-wound skin.



HIGH
ABSORBENCY¹

HIGH WET
STRENGTH

INCREASED
VERTICAL WICKING¹

SIZES AND CODES

ACTIVHEAL® AQUAFIBER™ EXTRA					
SIZE (CM)	QTY PER CARTON	QTY PER CASE	PRODUCT CODE	NHS SUPPLY CHAIN	DT PIP CODE
5x5	10	10	9023625	ELY795	407-5339
10x10	10	10	9023632	ELY796	407-5347
15x15	5	10	9023649	ELY797	407-5628
20x30	5	10	9023656	ELY798	407-5321
2.5x30.5 Ribbon	5	10	9023663	ELY799	407-5602
2x46 Ribbon	5	10	9023670	ELY800	407-5610

ACTIVHEAL® AQUAFIBER™ Ag INDICATED FOR THE MANAGEMENT OF INFECTED WOUNDS OR WOUNDS THAT ARE AT RISK OF INFECTION.

KEY FEATURES:

// High silver content // Fast initial elution of silver ions // Sustained silver release // High wet strength //



KEY BENEFITS

ActivHeal® AquaFiber™ Ag Antimicrobial wound dressing is a sterile non-woven pad consisting of a high M (mannuronic acid) calcium alginate and carboxymethylcellulose (CMC). Silver ions are released in the presence of wound fluid.


As fluid is absorbed, the alginate forms a soft gel which assists in maintaining a moist environment for optimum wound healing and allows intact removal. Silver ions released in the presence of wound exudate are an effective antimicrobial agent for up to 7 days, based on in-vitro testing², against a broad spectrum of microorganisms frequently associated with bacterial colonisation and infection of wounds.

INDICATIONS FOR USE

ActivHeal® AquaFiber™ Ag Antimicrobial wound dressing is indicated for the following moderately to heavily exuding, partial to full thickness wounds:

- ✓ **Post operative wounds**
- ✓ **Trauma wounds** (dermal lesions, trauma injuries or incisions)
- ✓ **Leg ulcers**
- ✓ **Pressure ulcers**
- ✓ **Diabetic ulcers**
- ✓ **Partial and superficial burns**
- ✓ **Graft and donor sites**
- ✓ **Cavity wounds**

SILVER SALTS ARE ATTACHED TO THE FIBERS AND IONIC SILVER IS RELEASED IN THE PRESENCE OF WOUND EXUDATE

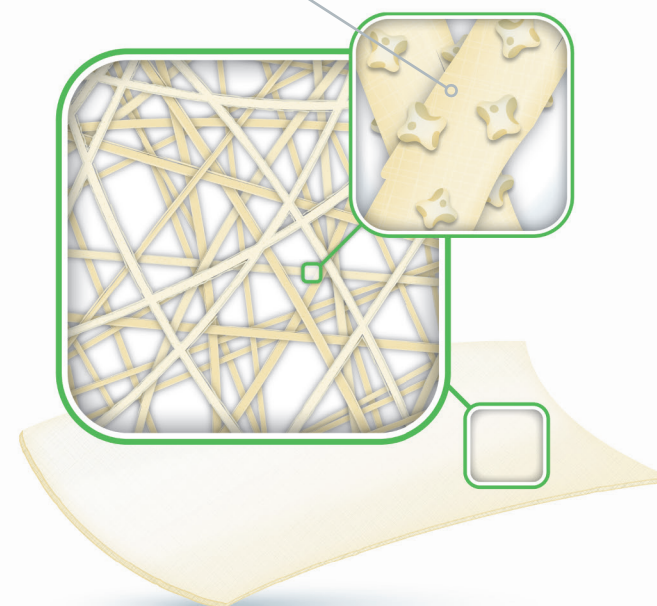

BACTERIA	AQUAFIBER™ AG			AQUACEL® AG		
	7 DAYS	14 DAYS	21 DAYS	7 DAYS	14 DAYS	21 DAYS
MRSE	✓	✓	✓	✓	✓	✗
MRSA	✓	✓	✓	✓	✓	✗
VRE	✓	✓	✓	✓	✓	✗
S.pyogenes	✓	✓	✓	✓	✓	✗
S.epidermidis	✓	✓	✓	✓	✓	✗
E.coli	✓	✓	✓	✓	✓	✓
C.albicans	✓	✓	✓	✓	✓	✗
P.aeruginosa	✓	✓	✓	✓	✓	✓

Fast initial and sustained anti-microbial activity²



AquaFiber™ Ag
16.43
g/cm²



AQUACEL® Ag
15.26
g/cm²

Absorbency¹
Higher absorbency can result in reduced dressing change frequency, minimising clinician time and disruption to the patient.



**HIGH SILVER
CONTENT**



**FAST INITIAL
ELUTION**



**SUSTAINED
SILVER RELEASE**



**HIGH WET
STRENGTH**



SIZES AND CODES

ACTIVHEAL® AQUAFIBER™ Ag				
SIZE (CM)	QTY	PRODUCT CODE	NHS SUPPLY CHAIN	DT PIP CODE
5x5	10	9023236	ELY778	405-6172
10x10	10	9023243	ELY779	405-6180
15x15	5	9023267	ELY780	405-6206
2.7x32 Ribbon	5	9023281	ELY781	405-6222



ACTIVHEAL® SILICONE FOAM DRESSINGS HAVE A GENTLE BUT SECURE ADHESION WITH MINIMAL EPIDERMAL STRIPPING OR PAIN ON REMOVAL.

KEY FEATURES:

// High total fluid handling // Perforated wound contact layer // Waterproof and bacterial barrier // Minimises epidermal stripping // Prevents pain at dressing changes // Protects the periwound

KEY BENEFITS

ActivHeal® Silicone Border and Non-Border Foam dressings are for use when adherence to the wound is a potential problem due to fragile skin. The perforated silicone coated wound contact later prevents the dressing from sticking to the wound and reduces pain during dressing changes. The perforations allow the uptake of exudate and prevents excess fluid causing maceration to healthy skin.



INDICATIONS FOR USE

For use on moderate to heavily exuding chronic and acute wounds. Suitable to use under compression bandaging, the dressing may be used throughout the healing process on:

- ✓ **Pressure ulcers**
- ✓ **Diabetic ulcers**
- ✓ **Venous and arterial leg ulcers**
- ✓ **Post operative surgical wounds**
- ✓ **Cavity wounds** (as a secondary dressing)
- ✓ **Skin abrasions**
- ✓ **Superficial and partial thickness burns**

SIZES AND CODES

ACTIVHEAL®		SILICONE FOAM BORDER			SILICONE FOAM NON-BORDER		
SIZE (CM)	QTY	PRODUCT CODE	DT PIP CODE	NHS SUPPLY CHAIN	PRODUCT CODE	DT PIP CODE	NHS SUPPLY CHAIN
5x5	10				9001548	402-5383	ELA847
7.5x7.5	10	9001470	402-5318	ELA843	9001555	402-5391	ELA834
10x10	10	9001487	402-5326	ELA838	9001562	402-5417	ELA844
12.5x12.5	10	9001494	402-5334	ELA840			
10x20	10	9001500	402-5342	ELA839	9001579	402-5441	ELA833
15x15	10	9001517	402-5359	ELA841	9001449	402-5425	ELA845
20x20	10	9001524	402-5367	ELA842	9001586	402-5433	ELA846
Sacral	10	9001531	402-5375	ELA848			

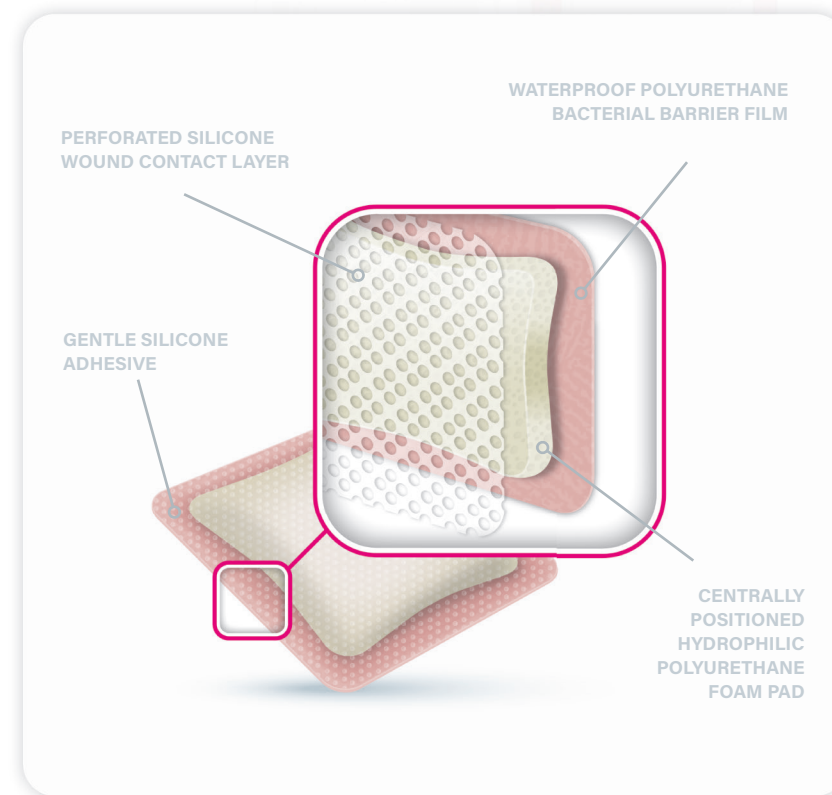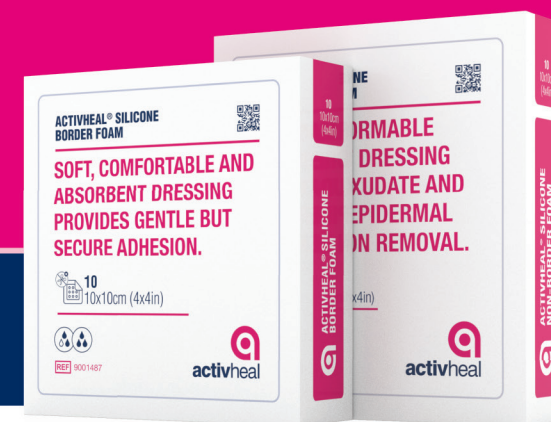

ACTIVHEAL® SILICONE FOAM LITE. THIN, SOFT AND DISCREET HYDROPHILIC FOAM PAD FOR NON AND LOW EXUDING WOUNDS.

KEY FEATURES:

// The dressing is easy to remove // Minimal epidermal stripping and pain during removal

KEY BENEFITS

ActivHeal® Silicone Foam Lite dressings helps to minimise maceration and damage to the surrounding skin. The perforated silicone wound contact layer aids atraumatic removal, enables exudate uptake, and prevents ingress of granulation tissue into the dressing. The dressing can be repositioned during application or lifted during wear time for observation.



INDICATIONS FOR USE

ActivHeal® Silicone Foam Lite dressing is indicated for use in the management of non to low exuding, chronic and acute wounds throughout the healing process. The dressing is indicated for use on the following wounds:

- ✓ **Pressure ulcers**
- ✓ **Diabetic ulcers**
- ✓ **Leg & foot ulcers**
- ✓ **Post-operative surgical wounds**
- ✓ **Trauma wounds** (including lacerations, abrasions, skin tears and blisters)
- ✓ **Superficial and partial thickness burns**
- ✓ **Radiation damaged skin**

SIZES AND CODES

ACTIVHEAL®		SILICONE FOAM LITE BORDER			SILICONE FOAM LITE NON-BORDER		
SIZE (CM)	QTY	PRODUCT CODE	NHS SUPPLY CHAIN	DT PIP CODE	PRODUCT CODE	NHS SUPPLY CHAIN	DT PIP CODE
5x5	10	9023687	ELA1135	412-9938			
5.5x12	10	9023694	ELA1136	413-0340			
7.5x7.5	10	9023700	ELA1137	413-0357	9023755	ELA1142	413-0423
10x10	10	9023717	ELA1138	413-0365	9023762	ELA1143	413-0431
15x15	10	9023724	ELA1139	413-0373	9023779	ELA1144	413-0415
10x20	10	9023731	ELA1140	413-0381	9023786	ELA1145	413-0407
10x30	10	9023748	ELA1141	413-0399			

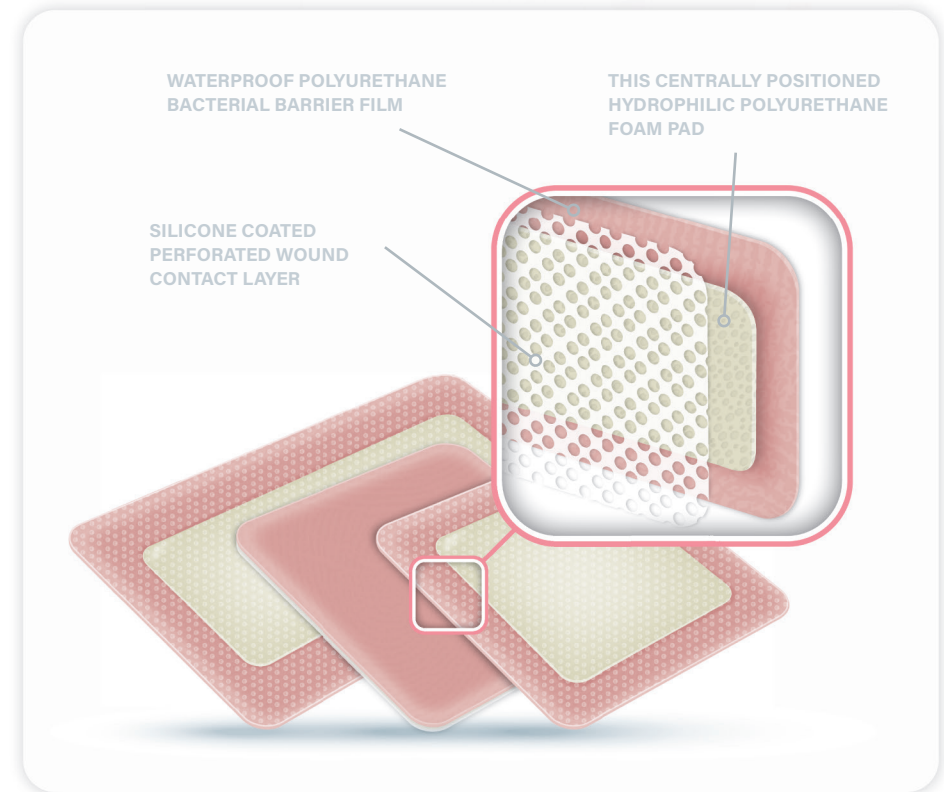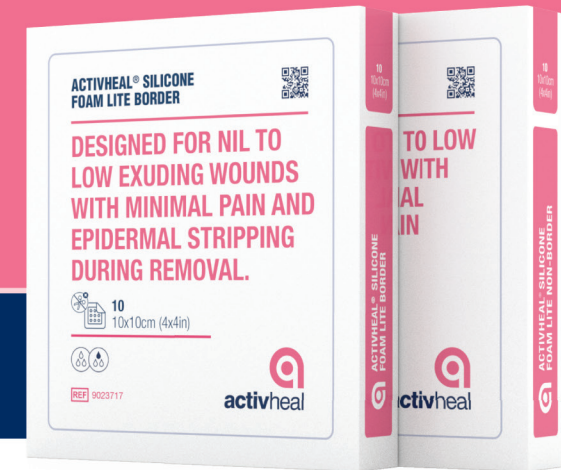

ACTIVHEAL® AQUAFIBER™ IS A SOFT, CONFORMABLE, HIGHLY ABSORBENT DRESSING THAT CONVERTS TO A SOFT CLEAR GEL WHEN IN CONTACT WITH WOUND EXUDATE

KEY FEATURES:

// High absorbency // Excellent wet and dry strength // Reduced lateral wicking // Clear gelling // Easy to use



KEY BENEFITS

ActivHeal® AquaFiber™ gels on contact with wound exudate and gently conforms to the wound surface. The reinforcement layer gives the dressing improved wet integrity providing the clinician with the confidence that the dressing can be removed intact. Calcium ions present within the alginate fiber of ActivHeal® AquaFiber™ exchange with sodium ions found in the wound exudate. The donation of calcium benefits the wound by stimulating cytokines which are essential to the wound healing process¹ and are also known to assist with the control of minor bleeding². ActivHeal® AquaFiber™ gives all the benefits of an alginate, while ensuring a moist wound environment is maintained.

INDICATIONS FOR USE

For use on moderately to heavily exuding chronic and acute wounds, including:


- ✓ **Pressure ulcers**
- ✓ **Leg ulcers**
- ✓ **Diabetic ulcers**
- ✓ **Cavity wounds**
- ✓ **Post surgical wounds**
- ✓ **Superficial and partial thickness burns**

SIZES AND CODES

ACTIVHEAL® AQUAFIBER™				
SIZE (CM)	QTY	PRODUCT CODE	NHS SUPPLY CHAIN	DT PIP CODE
5x5	10	10009145	ELY202	333-0701
10x10	10	10009146	ELY203	333-0719
15x15	5	10009147	ELY204	333-0743
2x46	5	10009148	ELY205	333-0750

NEEDED FIBER ENTANGLEMENT TO ENHANCE DRESSING INTEGRITY

HIDDEN INTEGRITY LAYER

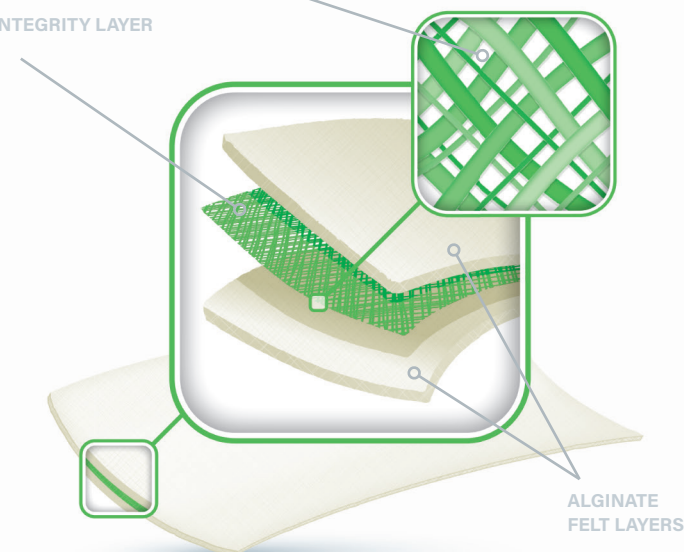

ACTIVHEAL® ALGINATE IS A PRIMARY SODIUM CALCIUM DRESSING THAT IS USED ON MODERATE TO HEAVILY EXUDING WOUNDS.

KEY FEATURES:

- // Excellent absorption of exudate // Reduces the risk of maceration // Aids autolytic debridement
- // Encourages granulation within the wound // Controls minor bleeding



KEY BENEFITS

ActivHeal® Alginate is manufactured by processing natural elements found in seaweed to produce felt and rope dressings. The absorbent alginate fiber in the dressing gel on contact with the wound fluid and gently conform to the wound surface.

Alginate dressings are known to donate calcium to the wound, assisting with the control of minor bleeding¹.

ActivHeal® Alginate is available in a diverse range of formats, the dressings are trusted by clinicians for the management of both topical and cavity wounds.

INDICATIONS FOR USE


For use on moderately to heavily exuding chronic and acute wounds, including:

- ✓ **Pressure ulcers**
- ✓ **Leg ulcers**
- ✓ **Diabetic ulcers**
- ✓ **Cavity wounds**
- ✓ **Donor sites**
- ✓ **Post surgical wounds**
- ✓ **Superficial and partial thickness burns**

SIZES AND CODES

ACTIVHEAL® ALGINATE				
SIZE (CM)	QTY	PRODUCT CODE	NHS SUPPLY CHAIN	DT PIP CODE
5x5	10	10007432	ELS139	301-6722
10x10	10	10007431	ELS140	301-6771
10x20	5	10007430	ELS141	301-6789
2x30 Rope	5	10007428	ELS142	301-6797


NEEDED FIBER ENTANGLEMENT TO ENHANCE DRESSING INTEGRITY

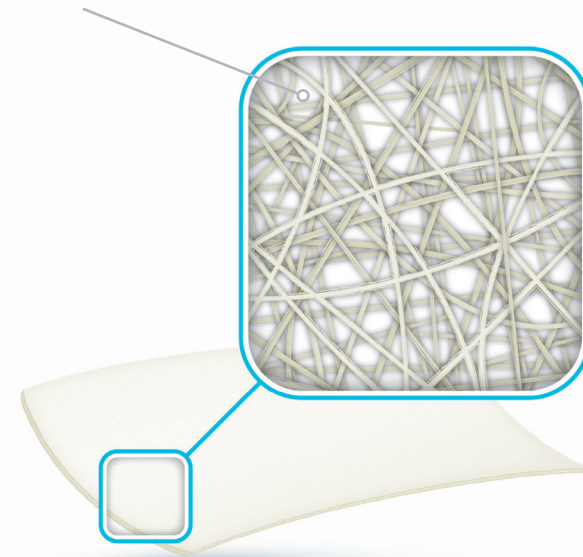

ACTIVHEAL® FOAM ADHESIVE IS A FLEXIBLE ABSORBENT FOAM DRESSING IDEAL FOR MODERATE TO HIGH EXUDATING WOUNDS.

KEY FEATURES:

// Effective total fluid handling // Low friction membrane // Waterproof and bacterial barrier //
Reduces the risk of maceration // Soft and comfortable

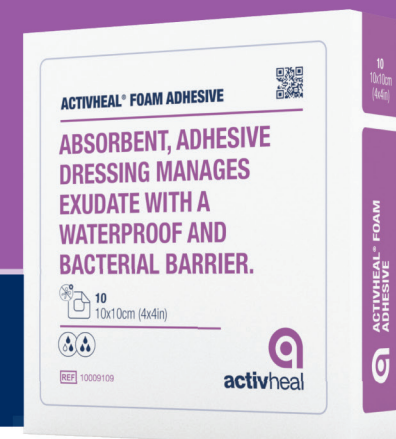

KEY BENEFITS

ActivHeal® Foam Adhesive provides secure adhesion to ensure the dressing remains in place allowing the patient to continue everyday activities confidently, whilst reducing the risk of leakage. The dressing conforms to the contours of the body which reduces the risk of rucking or catching on clothing and bedding.

INDICATIONS FOR USE

For use on moderate to heavily exuding chronic and acute wounds. Suitable for use under compression bandaging.


- ✓ Pressure ulcers
- ✓ Venous leg ulcers
- ✓ Diabetic ulcers
- ✓ Lacerations
- ✓ Abrasions
- ✓ Post operative surgical wounds
- ✓ Superficial and partial thickness burns
- ✓ Cavity wounds (as a secondary dressing)

SIZES AND CODES

ACTIVHEAL® FOAM ADHESIVE				
SIZE (CM)	QTY	PRODUCT CODE	NHS SUPPLY CHAIN	DT PIP CODE
7.5x7.5	10	10010232	ELM161	358-3390
10x10	10	10009109	ELA210	304-3767
12.5x12.5	10	10009110	ELA211	317-0495
15x15	10	10009111	ELA212	304-3775
20x20	10	10009112	ELA213	314-6065

WATERPROOF POLYURETHANE BACTERIAL BARRIER MEMBRANE


CENTRALLY POSITIONED HYDROPHILIC POLYURETHANE FOAM PAD

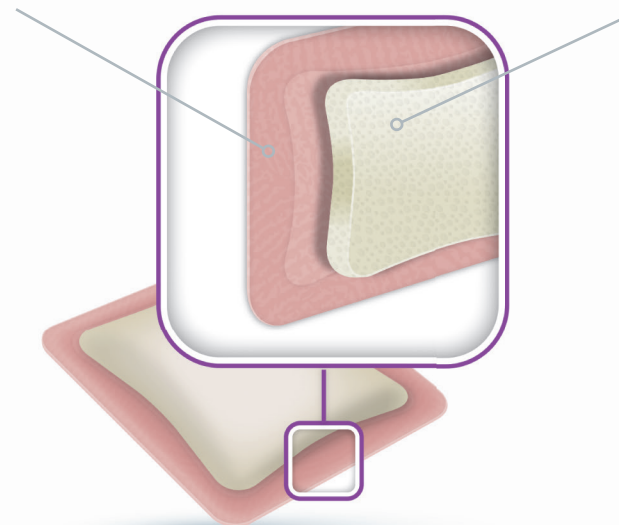

ACTIVHEAL® FOAM NON-ADHESIVE IS A VERSATILE DRESSING DEVELOPED TO OFFER PROTECTION AND ABSORBENCY ON MODERATE TO HEAVILY EXUDING WOUNDS.

KEY FEATURES:

// Effective total fluid handling // Low friction membrane // Bacterial barrier //
Soft and comfortable // Ideal under compression

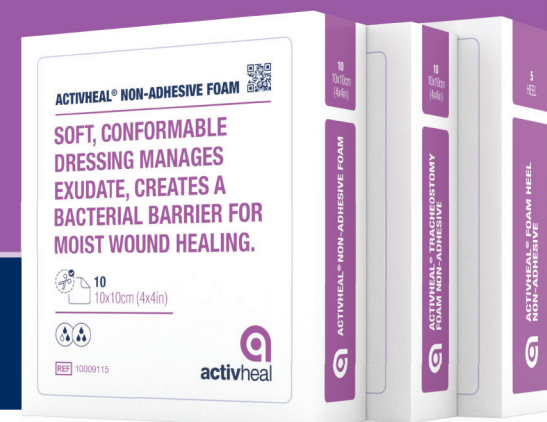

KEY BENEFITS

ActivHeal® Foam Non-Adhesive soft and flexible properties mean patients can wear it underneath clothing with comfort and confidence. The ActivHeal® Foam Tracheostomy dressing features a versatile fenestration, allowing it to fit neatly around a tracheostomy tube, intubation tube, cannula insertion, external bone fixators and appropriately sized wound drainage tubes. The ActivHeal® Foam Heel dressing, a two layer, soft pre-moulded dressing, specifically designed to apply to the heel and other bony anatomical areas.


INDICATIONS FOR USE

For moderate to heavily exuding chronic and acute wounds. Suitable for use under compression bandaging.

- ✓ Pressure ulcers
- ✓ Venous leg ulcers
- ✓ Diabetic ulcers
- ✓ Lacerations
- ✓ Abrasions
- ✓ Post operative surgical wounds
- ✓ Superficial and partial thickness burns
- ✓ Cavity wounds (as a secondary dressing)

SIZES AND CODES

ACTIVHEAL® FOAM NON ADHESIVE				
SIZE (CM)	QTY	PRODUCT CODE	NHS SUPPLY CHAIN	DT PIP CODE
5x5	10	10009113	ELA214	314-6073
7.5x7.5	10	10009114	ELA245	
10x10	10	10009115	ELA216	304-3791
10x20	10	10011265	ELA246	317-0503
20x20	10	10009117	ELA218	314-6081
18x12 Heel	5	10010227	ELM160	356-1131
10x10 Tracheostomy	10	10009118	ELM335	382-2483

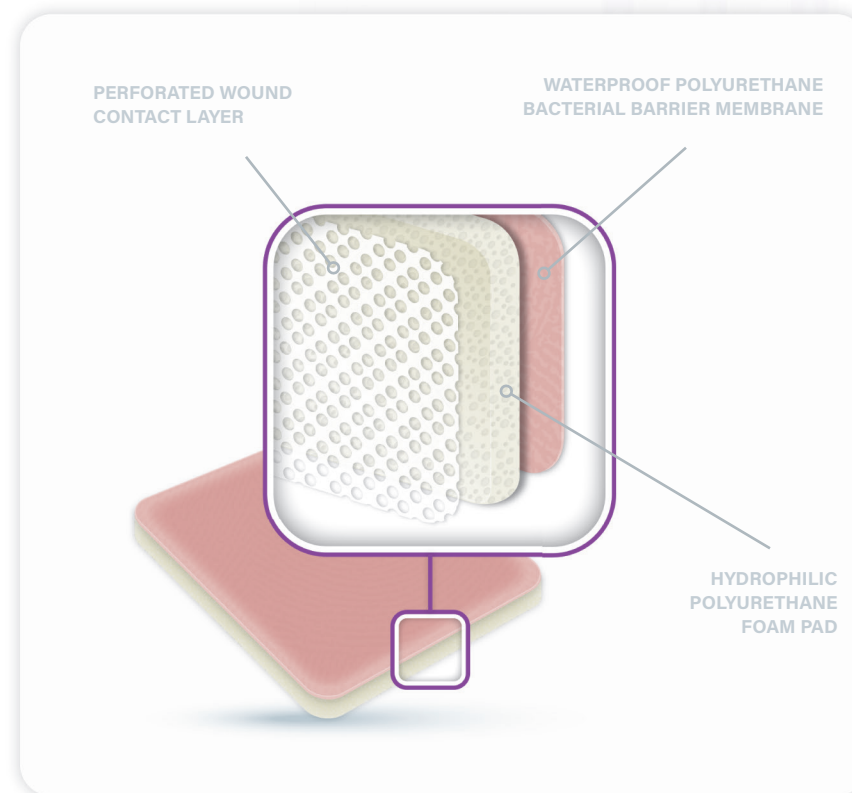

ACTIVHEAL® SILICONE WOUND CONTACT LAYER MINIMISES TRAUMA TO THE WOUND BED AND PERIWOUND SKIN.

KEY FEATURES:

// Soft silicone coating // Can remain in place for up to 14 days // Maintains a moist wound environment // Highly conformable // Minimises epidermal stripping // Minimises pain at dressing changes // Protects the periwound



KEY BENEFITS

ActivHeal® Silicone Wound Contact Layer dressing is a primary dressing for use when adherence to the wound is a potential problem. The dressing is a soft silicone wound contact layer with a silicone coating on both sides allowing exudate to pass through into the secondary dressing. ActivHeal® Silicone Wound Contact Layer protects the sensitive granulated wound bed and ensures a moist wound environment. It can remain in place for up to 14 days, during which the secondary outer absorbent pad can be changed without damaging the primary silicone wound contact layer.

INDICATIONS FOR USE


Nil to heavily exuding chronic and acute wounds. Can be used in conjunction with Negative Pressure Wound Therapy (NPWT). The dressing may be used throughout the healing process on:

- ✓ Skin tears
- ✓ Surgical incisions
- ✓ Second degree burns
- ✓ Pressure ulcers
- ✓ Skin grafts and donor sites
- ✓ Venous and arterial ulcers

SIZES AND CODES

ACTIVHEAL® SILICONE WOUND CONTACT LAYER				
SIZE (CM)	QTY	PRODUCT CODE	NHS SUPPLY CHAIN	DT PIP CODE
5x7.5	10	9001241	ELA849	399-7459
10x10	10	9001258	ELA835	399-7442
10x20	10	9001265	ELA836	399-7467
15x15	10	9001272	ELA837	399-7475

PERFORATED HOLES TO ALLOW
EXUDATE TO PASS THROUGH TO A
SECONDARY DRESSING

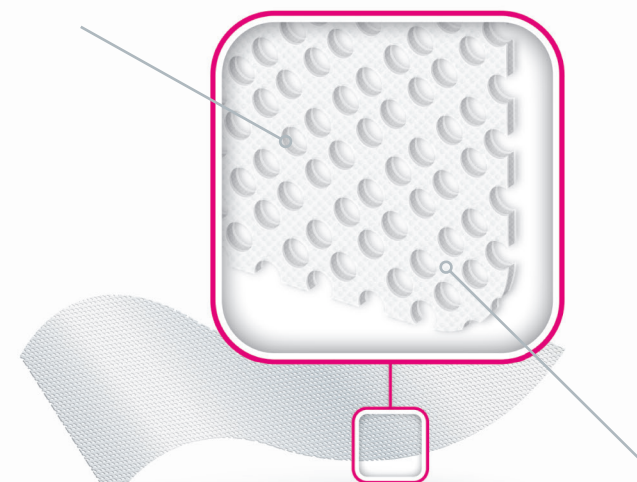

KNITTED FABRIC COATED BOTH
SIDES WITH SOFT SILICONE

ACTIVHEAL® HYDROGEL SPECIFICALLY FORMULATED WITH HIGH WATER CONTENT, THIS HYDROGEL DONATES MOISTURE TO DRY, NECROTIC AND SLOUGHY WOUNDS.

KEY FEATURES:

- // High fluid donation // Aids autolytic debridement // Rehydrates dry necrotic wounds
- // Removes devitalised tissue



KEY BENEFITS

ActivHeal® Hydrogel is a crucial product in any wound care treatment when fluid donation is essential to deliver important rehydration to a wound. ActivHeal® Hydrogel encourages healing by providing a moist wound healing environment thereby supporting autolytic debridement and desloughing. Highly viscous ActivHeal® Hydrogel conforms to any wound surface due to its amorphous state.


INDICATIONS FOR USE


ActivHeal® Hydrogel is used as a primary dressing indicated for use on dry and sloughy wounds with nil to low exudate:

- ✓ Pressure ulcers
 - ✓ Leg ulcers
 - ✓ Diabetic ulcers
- ✓ Lacerations and abrasions
 - ✓ Graft and Donor sites
 - ✓ Cavity wounds

SIZES AND CODES

ACTIVHEAL® HYDROGEL				
SIZE	QTY	PRODUCT CODE	NHS SUPPLY CHAIN	DT PIP CODE
8g	10	10015047	ELA639	377-9956
15g	10	10015048	ELG018	301-6847

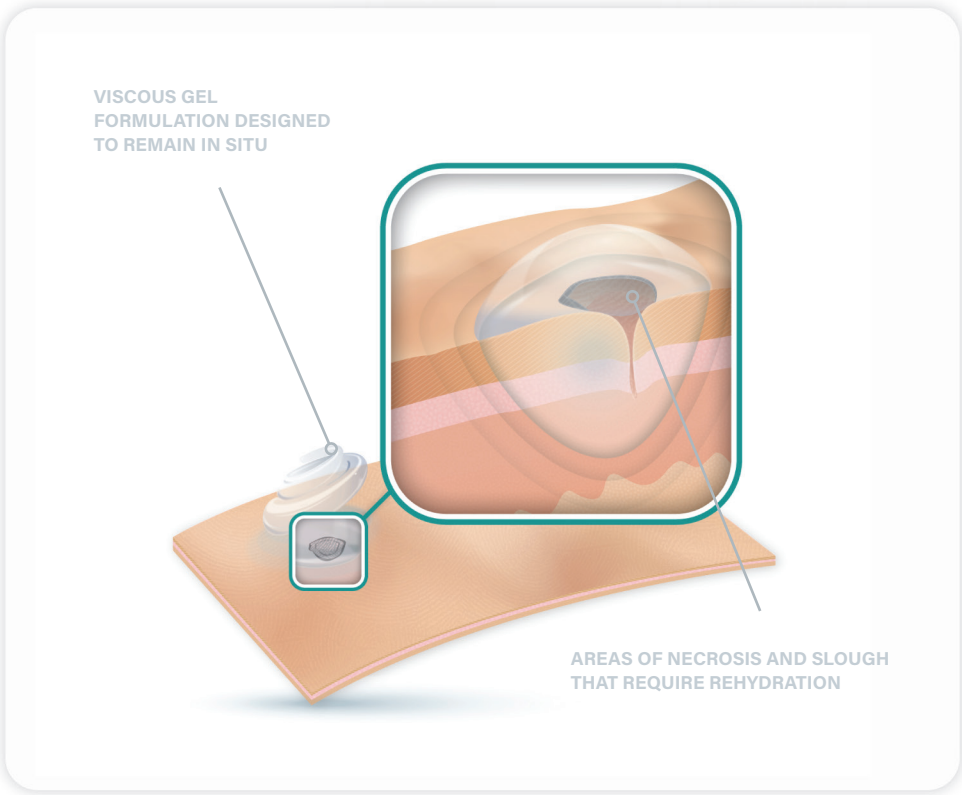

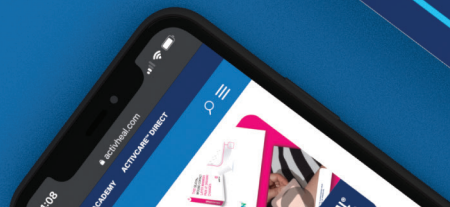



activheal.com



Advanced Medical Solutions Ltd

Premier Park, 33 Road One, Winsford Industrial Estate, Winsford, Cheshire, CW7 3RT, United Kingdom

T 08444 125 755 E activheal@admedsol.com W activheal.com

ActivHeal its logo, AquaFiber and the Advanced Medical Solutions logos are trade marks of Advanced Medical Solutions Ltd. Copyright © 2021.

